



CID-LX



Reliable
Simple
Standalone

MADE IN FRANCE

PORTABLE
RESPIRATORY
POLYGRAPH

www.cidelec.net





Benefit from the
combination of
CIDELEC technology
and the precision of
PneaVoX



Signal readings and
analyses

**OUR OBJECTIVE : ENABLING YOU TO USE ACCURATE,
COMPLETE AND RELIABLE RECORDINGS.**

STANDALONE

PORTABLE

SIMPLE

RELIABLE

COMFORTABLE

/ HOW IT WORKS



CID-LX

- **Installed** on the patient
- **12 channels** integrated into the device
- **Use** in polygraph




CIDELEC software

- **Automatic analysis** of the signals
- **Review** of the traces
- **Archiving** of data
- Customisable **summaries**

CLICK 'N CID



**DOWNLOAD THE APP
TO SEE A 3D INSTALLATION
OF OUR PRODUCTS.**

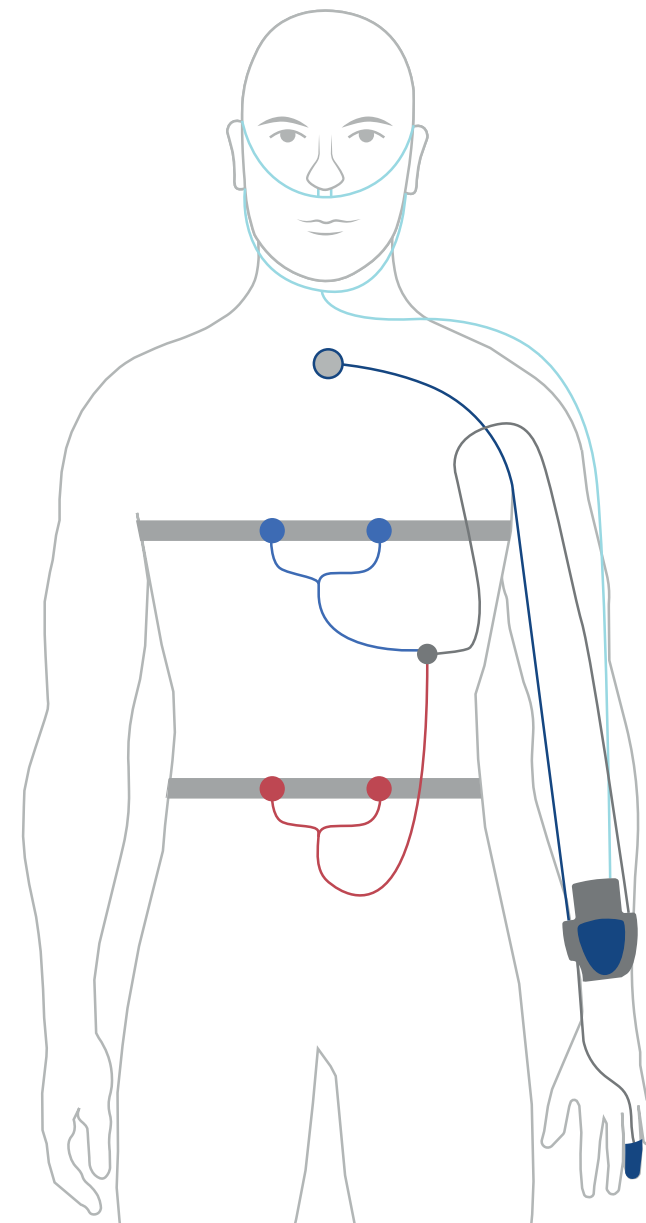

/ WE OFFER

CIDELEC provides:

- **Training and installation of products on site** by our team,
- **Innovative technology of our systems** to obtain precise and reliable measurements
- **Analysis and processing of information** with the creation of personalised reports,
- **After-sales service**, technical assistance and technical expertise



**DEVICES DESIGNED
& MADE IN FRANCE**

**CIDELEC,
30 YEARS
YOUR PARTNER**

/ PERFORMANCE & QUALITIES

For almost 30 years, our devices to aid the diagnosis of sleep-related or sleep-aggravated pathologies have been designed and manufactured in France.

CIDELEC supports you throughout their use: presentation, sales, installation, user training, telephone assistance, after-sales service.

The CID-LX-206d model also has a pressure channel for the connection of a pneumotachograph when the patient uses cPAP/BIPAP treatment.

12 channels available

Channel 13: only on CID-LX-206d

Dimensions: 32 x 82 x 114 mm- Weight: 135 g- Battery: Li-Po 1700 mAh - 3.7V

Technical characteristics CID-LX

CHANNELS	BANDWIDTH	SAMPLING FREQUENCY	STORAGE	PRECISION	POINTS	ELONGATION	OTHER
Breathing sound	200 - 2000 Hz	4000 Hz	Sound intensity at 16 Hz		256		Sensitivity 20 - 80 dB Adaptive threshold
Snoring	20 - 200 Hz	4000 Hz	Sound intensity at 16 Hz		256		Sensitivity 60 - 120 dB Threshold 76 dB
Suprasternal pressure	0.02 - 20 Hz	4000 Hz	8 Hz		4096	+/- 100 Pa	
Position		1 Hz	1 Hz				5 positions
Actimeter		1000 Hz	8 Hz				
Nasal flow	0 - 10 Hz	4000 Hz	256 Hz		65536	+/- 300 Pa	
Machine pressure	0 - 10 Hz	4000 Hz	256 Hz	+/- 25 Pa	4096	0 - 2 kPa	Up to 4 kPa on request
SpO₂⁽¹⁾			75 Hz	+/- 3% 100% ⁽²⁾	100	0 - 100%	Averaged over 4 pulse cycles
Pulse rate⁽¹⁾			8 Hz	+/- 5 BPM ⁽²⁾	128	40 - 240 BPM ⁽²⁾	
Photoplethysmogram⁽¹⁾			75 Hz				
Inductive straps	0.1 - 10 Hz		32 Hz		65536		
Pneumotachograph⁽³⁾	0 - 10 Hz	4000 Hz	16 Hz	+/- 4%	4096	+/- 1 litre/s	

(1) NONIN manufacturer

(2) Under the least favourable conditions

(3) Only available on CID-LX-206d

/ PNEAVOX

PneaVoX technology is unique.

One sensor records 3 physiological parameters :

- Buccal and nasal **breathing**,
- **Respiratory effort** via suprasternal pressure to differentiate between obstructive, central and combined apneas,
- **Snoring** (energy, intensity).

Finally, the PneaVoX sound sensor **analyses upper airway resistance** by measuring the sound intensity.

"The PneaVoX sound sensor, to improve differentiation between sleep disorders via the analysis of tracheal sounds."

PNEAVOX[®]
TECHNOLOGY

/ SCIENTIFIC BIBLIOGRAPHY

M. Glos, A. Sabil, K.S. Jelavic, C; Schöbel, I. Fietze, T. Penzel. *Characterization of respiratory events in obstructive sleep apnea using suprasternal pressure monitoring.* *J Clin Sleep Med.* 2018; 14(3): 359-369.

Penzel T, Sabil A. *The use of tracheal sounds for the diagnosis of sleep apnoea.* *Breathe* 2017; 13: e37-e45.

A. Amaddeo, M. Fernandez-Bolanos, J.O. Arroyo, S. Khirani, G. Baffet, B. Fauroux. *Validation of a Suprasternal Pressure Sensor for Sleep Apnea Classification in Children,* *Journal of Clinical Sleep Medicine, Vol. 12, No. 12, 2016.*

[...]

PURCHASE



CIDELEC

20 rue des Métiers

49130 SAINTES GEMMES SUR LOIRE - FRANCE

Tel: +33 (0)2 41 66 20 88 - Fax: + 33 (0)2 41 79 07 76

Email: service.commercial@cidelec.net

Website: www.cidelec.net

The content of this brochure (text, illustrations, drawings, images, etc.) is the exclusive property of CIDELEC. The CIDELEC logo and names, CID, CID102, TrackBox and PneaVoX are trademarks registered by CIDELEC. CIDELEC reserves the right to modify this document at any time.

The CID-LX is a class IIa medical device, manufactured by CIDELEC - CE No. 0459

The CID-LX is a device for collecting physiological signals for the diagnosis of sleep disorders.

Read the product instructions carefully before use. Document modified on 01/2022